# INOD

## Q.P. FIBER POST PLUS


QUALIFIED PRODUCT

Step 1

Step 2

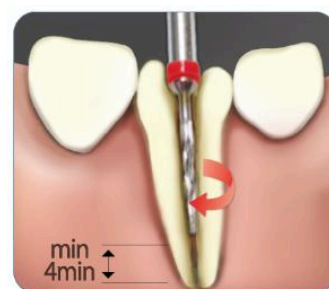
Step 3


1. Root canal filled with Gutta-percha.
2. Remove root filling using selected universal drill(1000~2000rpm).  
Leave minimum 4mm of filling apically.
3. Select post size and matching drill. Use appropriate drill to widen and shape the root canal.

Step 4

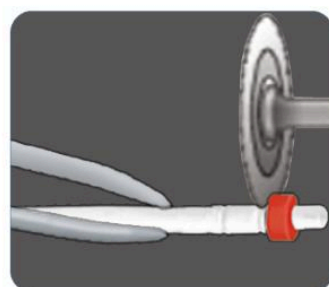
Step 5

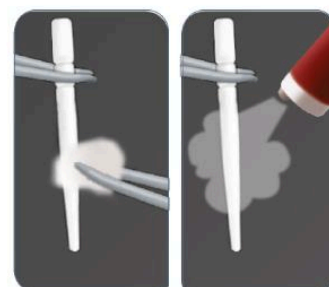Step 6

4. Insert post to check the fit in the prepared canal. Should be able to easily remove the post.
5. Shorten post to the required length using a diamond disk or diamond bur.
- \* Note) Do not use scissors, scalpels or similar tools since they will crush the post and ruin the structure and stability.
6. Disinfect the post with alcohol and dry it with water-free and oil-free air.

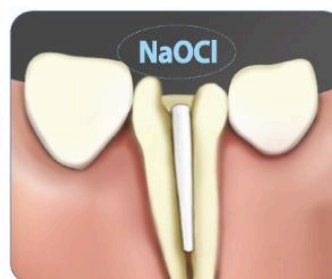Step 7


Step 8

7. Clean or disinfect the root canal with 2.5-5.25% sodium hypochlorite solution(NaOCl). Rinse immediately with water and dry with paper points.

\* Note) Do not use hydrogen peroxide(H<sub>2</sub>O<sub>2</sub>), other disinfectants, and EDTA solutions because their residues can impair the adhesion strength and curing reaction of cement.

8. Begin cement dispensing in the apical root canal. Keep elongation tip immersed in the cement initially and during the entire cement application process to avoid trapped air in the cement.

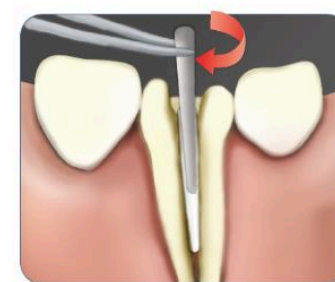Step 9

Step 10

Step 11

9. Seat the post immediately. Twist slightly and apply moderate pressure to hold in position.
10. Apply moderate pressure while removing excess with appropriate instruments or a cotton pellet.
11. Light cure cement or allow self-cure as cement manufacturer recommendation.

Model	Diam(Ø, mm)	Point Diam(Ø, mm)	Length(Ø, mm)	Drill Diam(Ø, mm)
QPA14	1.14	0.6	20	1.14
QPA25	1.25	0.7	20	1.25
QPA55	1.55	0.8	20	1.55
QPA9	0.9	0.7	20	0.9
QPA10	1.0	0.8	20	1.0
QPA20	1.2	0.65	20	1.2
QPA40	1.4	0.75	20	1.4

(mm)

### ► Side effect

There is no special side effect.

### ► Caution :

- Avoid exposure to direct light during use.
- Avoid getting dirty.

### ► Storage condition

Keep it at room temperature.

### ► Contraindication :

It does not intend to treat children or pregnant or nursing women.

- Make sure to use diamond disk to cut.
  - Use sponge with alcohol to clean
- Prevent direct contact with water and oil.

### ► Disposal

It should be handled according to local laws and regulations.

## INNO DENTAL CO., LTD.

812, Byucksan Digital Valley, 6-cha, 219,  
Gasandigital 1-ro, Geumcheon-gu, Seoul 08501, Republic of Korea  
Tel : +82-2-757-1112 / Fax : +82-2-757-9990 / www.inod.co.kr  
SRN Number : KR-MF-000041291

EC REP

Name : S.B. PHARMA GMBH  
Address : Bunsenstr, 14 53121 Bonn, Germany  
Tel : +49(0) 228 52266916

CE  
0068

# Instructions for Use

This product is a non-metallic root canal post as a fiber post. After curing root canal, insert it into the root of a tooth in order to be easy core forming. This post is intended to protect teeth thanks to distracting stress that can effect to residual teeth supporting while recovering of artificial root canal.

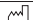
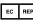
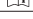
**- Product Name : Non-metallic root canal post**

**- Brand name : Q.P. Fiber Post Plus Refill**

## Indicate for use

Non-metallic filler used during root canal operation to connect the root and support it.

### Explanation of symbol

	This symbol provides notification of important information. When this symbol is on a device, it refers to important information in the accompanying documents.
	This symbol is Date of manufacture
	This symbol is manufacture
	This symbol is EU representative.
	This symbol is lot number
	This symbol is use by
	This symbol is do not re-use
	This symbol is "DO NOT USE IF PACKAGE IS DAMAGED".
	Contents
	This symbol is "Consult instruction for use".

### Product Description

Model	Diam(ø, mm)	Length(ø, mm)	Drill(ø, mm)
QPA9	0.9	20	0.9
QPA10	1.0	20	1.0
QPA14	1.14	20	1.14
QPA20	1.2	20	1.2
QPA25	1.25	20	1.25
QPA40	1.4	20	1.4
QPA50	1.5	20	1.5
QPA55	1.55	20	1.55
QPA80	1.8	20	1.8
QPA200	2.0	20	2.0

S.B. PHARMA GMBH  
Bunsenstr. 14 53121 Bonn, Germany  
TEL: +49 (0) 228 52266916

## How to use

- Step.1 Root Canal filled with Gutta-percha.
- Step.2 Remove root filling using selected universal drill(1000~2000rpm). Leave minimum 4mm of filling apically.
- Step.3 Select post size and matching drill. Use appropriate drill to widen and shape the root canal.
- Step.4 Insert post to check the fit in the prepared canal. Should be able to easily remove the post.
- Step.5 Shorten post to the required length using a diamond disk or diamond bur. Note: Do not use scissors, scalpels or similar tools since they will crush the post and ruin the structure.
- Step.6 Disinfect the post with alcohol and dry it with water-free and oil-free air.
- Step.7 Clear or disinfect the root canal with 2.5-5.25% sodium hypochlorite solution(NaOCl). Rinse immediately with water and dry with paper points. Note: do not use hydrogen peroxide (H2O2), other disinfectants, and EDTA solutions because their residues can impair the adhesion strength and curing reaction of cement.
- Step.8 Begin cement dispensing in the apical root canal. Keep elongation tip immersed in the cement initially and during the entire cement application process to avoid trapped air in the cement.
- Step.9 Seat the post immediately. Twist slightly and apply moderate pressure to hold in position.
- Step.10 Apply moderate pressure while removing excess with appropriate instruments or a cotton pellet.
- Step.11 Light cure cement or allow self-cure as cement manufacturer recommendation.

### Contraindication

It does not intend to treat children or pregnant or nursing women

### Storage condition

Keep it at room temperature.

### Disposal

It should be handled according to local laws and regulations.

### Side effect

There is no special side effect.

### Caution

1. Avoid exposure to direct light during use.
  2. Make sure to use diamond disk to cut.
  3. Avoid getting dirty.
  4. Use sponge with alcohol to clean.
- Prevent direct contact with water and oil.

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## INNODENTAL CO.,LTD